

Quick start guide

EN

Hermes medical engineering S.L.
www.minivap.com




1.

Connect the charger to the charging base (1.1) to charge the battery or adapter (1.2).


Attach (1.3) the MiniVAP to the battery or adapter.



1.2

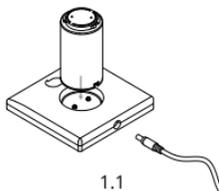


1.1

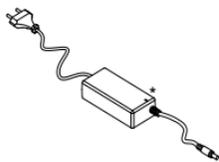

1.3


LED CHARGING LIGHTS*		
Color	Adapter	Battery
Green	Inactive	Charged
Red	Heating	Charging

When the MiniVAP is switched on, the front LED (temperature) (2.1) and black LED (battery level) (2.2) will both light up.

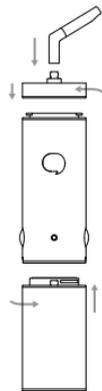


2.

Switch on the Vaporizer by pressing the two buttons for 2 seconds.

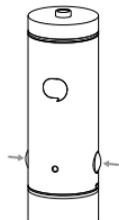

The front LED will flicker to indicate that the device is heating up.

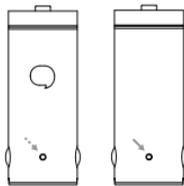


2.1

2.2

3.

Push the buttons '+' or '-' to regulate the temperature.

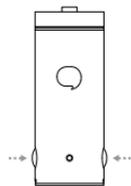

FRONT LIGHT	
Color	Temperature*
Blue	365 °F / 185 °C
Pink	392 °F / 200 °C
Green	419 °F / 215 °C
Red	446 °F / 230 °C

**aprox. value, (depends on the depth or duration of the inhalation).*


BACK LIGHT		
Color	Battery	Life*
Green	Complete	2-3 hours
Orange	Medium	1-2 hours
Red	Low	1 hour

**aprox. value, (depends on the number of inhalations and the model).*

4.

Fill the basket with the (previously finely shredded) plant or place the extract on the mesh filter (4.1).

Install the O-ring filter inside the flexicone+ or flexicone and then insert the basket or mesh filter inside it (4.2).

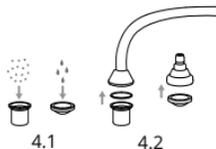


4.1

4.2


Introduce the flexicones through the opening of the lid. Place the lid on the MiniVAP and screw it up gently.

Use the MV mouthpiece for the flexicone+. Install the elbow in the mouthpiece for the flexicone.

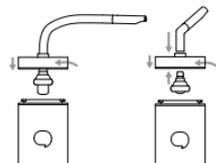

5.

Wait until the front light has stopped flickering.

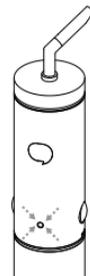

Then inhale deeply and continuously for at least 3 seconds.

The amount of vapor depends on the depth and duration of the inhalation.


Slowly exhale enjoying the taste and aroma of a high quality vaporizer.

You can repeat this step between 5 to 15 times depending on the concentration and amount of plant/extract.

6.

Turn off the vaporizer by pressing the buttons for 2 seconds.

When the vaporizer has cooled down, vaporized plant material can be cleaned from the loading chamber, basket and filter by using the cleaning brush.

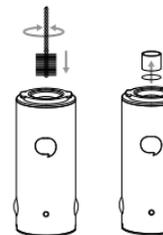

Useful Tips

The **basket** must not be filled to the brim, so that the hot air can freely circulate and it is easier to inhale.

Retain the vapor in the lungs for a few seconds for greater absorption of the active principle.

IMPORTANT

The front LED may flicker if the inhalation is deep or until it reaches the previously selected temperature.

MiniVAP saves energy. If no button is pressed or if not inhaled within 5 minutes, the vaporizer temperature falls to a minimum (blue).

If the buttons of the vaporizer are not pressed, it automatically switches itself off after 10 minutes.

For in-depth **cleaning** immerse the filter, basket, and mouthpieces in alcohol and then rinse thoroughly with water.

Safety measures

Use and store the device in a dry place and keep protected from the elements.

We recommend using the device in well ventilated areas for proper cooling. The device must not be covered with any materials that restrict direct contact with the external environment.

Only use the device to inhale and never exhale air into the vaporizer as this can cause damage due to overheating, deterioration of the battery, and caking and waste inside internal parts.

Place the device and the charging base on a horizontal and flat surface, away from any source of humidity, heat, and inflammable gases.

Keep the device in sight at all times and out of reach of children.

Always use the charger supplied by the manufacturer, when connecting it to the vaporizer and for recharging the battery.

If the vaporizer is damaged or malfunctioning, please contact the manufacturer.

support@minivap.com

IMPORTANT

Store the vaporizer in a safe place only when it is turned off.

There is a risk of burns in contact with the vaporization chamber.

Never introduce an object or liquids object or liquid into the vaporization chamber.

Ensure that the battery is on charge when the device is switched off.

Take special care when placing the battery on wet surfaces, as it might lead to a short circuit.

Fully charge and discharge the battery at least once a month.

Clean the inside of the vaporizer when it has cooled down and the unit disconnected from the battery or the charger.

For more information, check the Manual or FAQs at the Support menu of our website:

www.minivap.com

